




BODY HEALTH[®]
Shape your Future

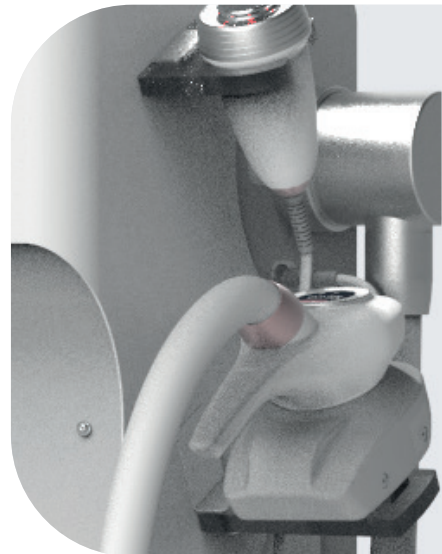

UNYQUE

Registration ANMAT: PM 1886-1
Model: BHS 202

INNOVATIVE AND EXCLUSIVE DESIGN FOR A BETTER USER EXPERIENCE

The new Body Health line with a movable head and 14" LED touchscreen makes treatment interaction and application easier.

In addition to its cutting-edge technology, UNYQUE ensures ultimate results.



TOUCHSCREEN

With pivot system.



APPLICATION SOFTWARE

Intuitive multilingual design.



WEB MANAGEMENT PLATFORM

24h remote geolocation and management.



TRACK SCREEN SYSTEM

Parameters are displayed directly on the handpiece.




SOFTWARE




Interactive Software




Gender Selection

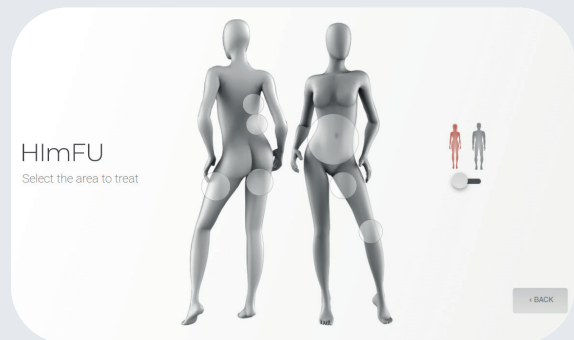

Alarm system for better accuracy in the treatment of quadrant

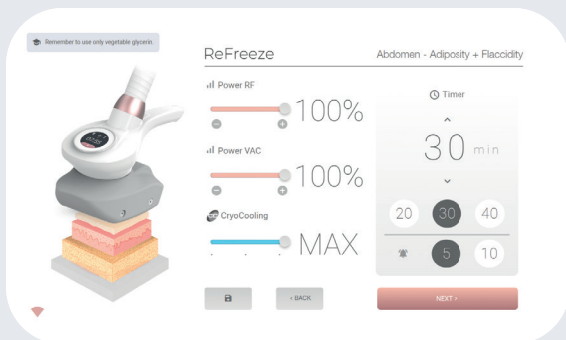

Selection and modification of parameters

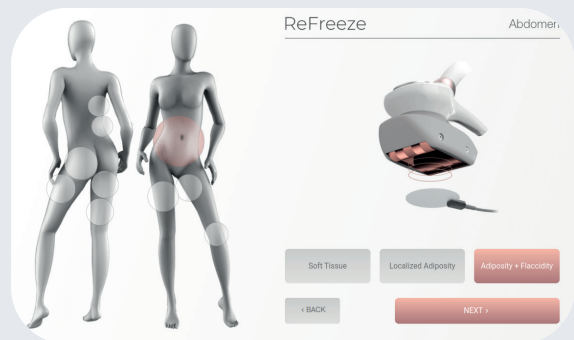

THERAPIES

The market's most complete equipment that applies therapies with TrackScreen technology and a new management software, delivering amazing results. In a single unit, this equipment combines ReFreeze, CryoFrequency and HImFU therapies.



ReFreeze ^{NUEVA}

Indicated for

- Lymphatic drainage.
- Cellulite.
- Tightening.
- Body contouring
- Localized adiposity



CryoFrequency

+ Multi and Monopolar
Radiofrequency for Face y Body

Indicated for

- Tightening.
- Cellulite.
- Body contouring.
- Localized adiposity.
- Wrinkles.



HImFU Corporal

Indicated for

- Localized adiposity.
- Body contouring.




ReFreeze NEW

· A therapy that combines two technologies in the same handpiece, **Unipolar + Multipolar Cryofrequency** together with simultaneous and continuous suction, which allow heat to penetrate faster and deeper into the treated area.

· Indicated to treat areas with generalized **adiposity, cellulite and flaccidity** while performing **lymphatic drainage**. Eliminates waste, toxins and fluid overload and achieves **instant and lasting effects**.

· In addition to **Body Contouring** for body shaping, the **CryoCooling** system easily applies cold on the skin, providing an anesthetic effect for greater patient comfort.

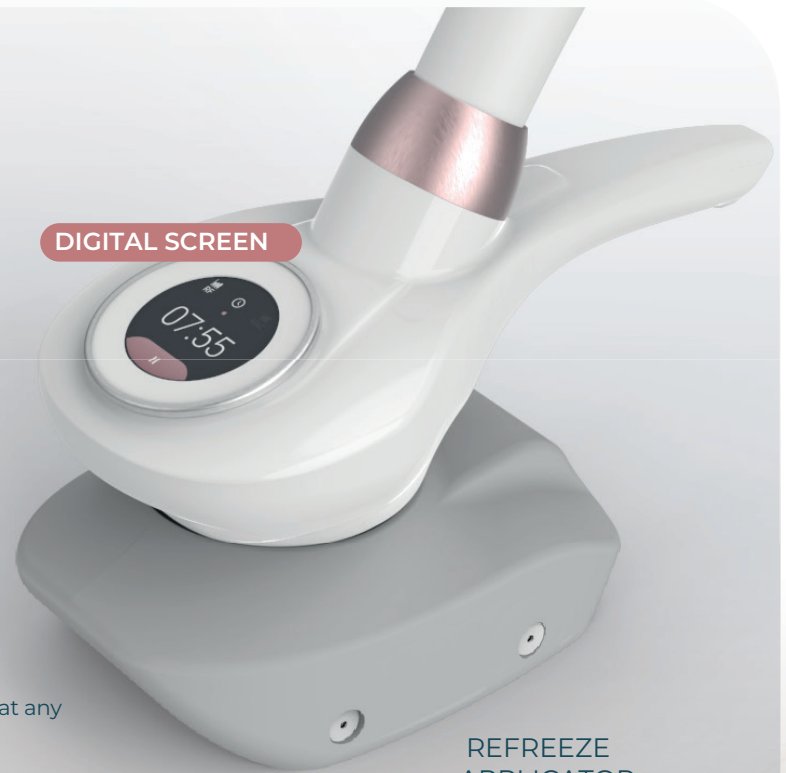

REFREEZE APPLICATOR




Benefits

- Progressive and long-lasting effect.
- May be applied on any skin type and at any time of the year.
- No hospitalization or recovery period required.
- Reduces volume and shapes the body.
- Painless sessions.
- Quick and easy to apply.
- No need for anesthesia




CryoFrequency

CryoFrequency is a **safe and non-invasive** therapy that combines the transmission of cold and radiofrequency, causing an internal heat in body tissues (over 55°C) that provides treatment for general signs of flaccidity and adiposity. The combination of **cold and heat** causes a **thermal shock** that stimulates the synthesis of new collagen and elastin fibers producing an **instant, progressive and long lasting lifting effect from the first session**. In addition, the unipolar mood generates an intensive heating of the fat cells in the hypodermis, activating the **apoptosis** of the adipocyte, thus reducing the adipose panniculus.

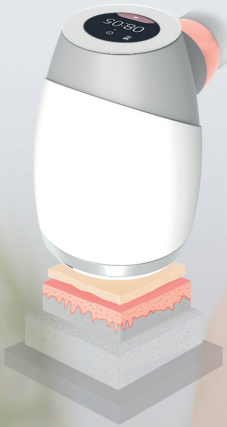

Benefits

- Progressive and lasting effect.
- May be applied on any skin type.
- Reduces volume and shapes the body.
- Painless sessions.
- Quick and easy to apply

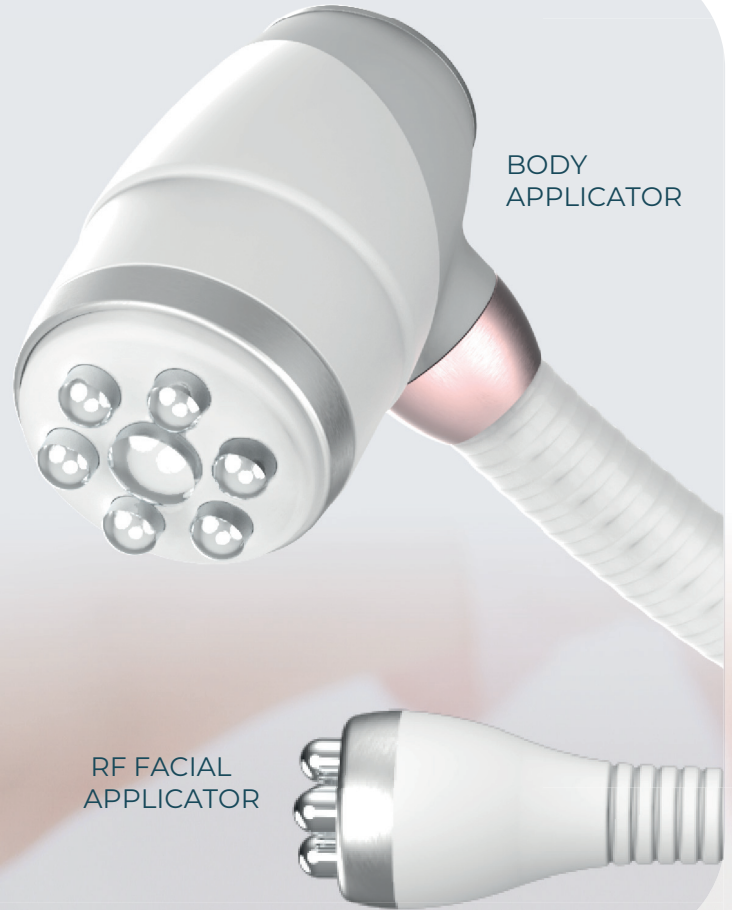


Multifrequencial Radio frequency

The Orbicular applicator device uses Tripolar technology with 4 poles that deliver a total of 400 watts generating crossed fields, which increases the temperature at the dermal level. This applicator can also be used with bipolar RF.

BODY HEALTH's CryoFrequency is unique. It has 6 surgical steel poles that generate three crossed power arcs, increasing the internal temperature with a controlled depth at the dermal level.

In addition, a lipolytic effect can be generated with the unipolar program. Since adipose tissue is resistant to power flow, the radiofrequency waves cause selective heating of this tissue, resulting in apoptosis of the adipose cells in the treated area.

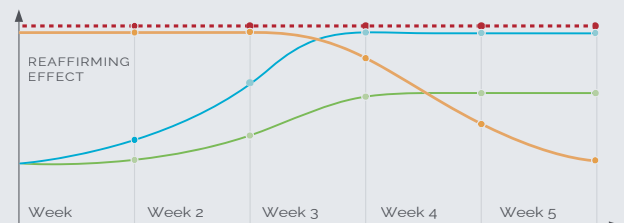


- Resultado instantáneo y duradero
- CryoFrequency Collagen Development
- Thermal Shock
- Other Radiofrequencies Collagen Development



HimFU Corporal

+ Dermoleds


HimFU Corporal is an advanced, innovative, unique therapy for the treatment of localized adiposity. Multifocal ultrasound waves are used; this results in a thermomechanical effect on adipose tissue causing its destruction. Patients can reduce between 1 and 4 cm of contour per treatment session obtaining an incredible body reshaping effect.

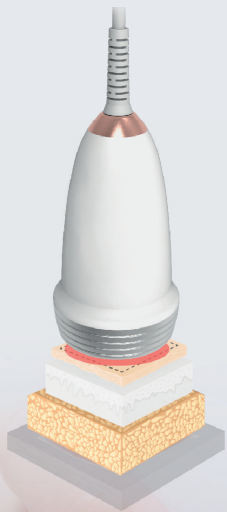

Benefits

- Not invasive.
- Painless sessions without the risk of burns.
- Safe and effective.

DERMOLEDS

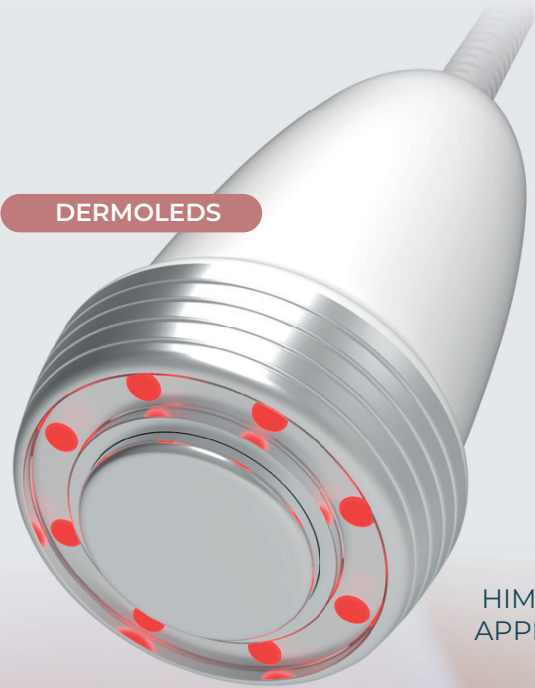

HIMFU
APPLICATOR



Applicator Frequency

The HimFU technology has a transducer that operates at 2.9 Mhz delivering MULTIFOCAL energy at a depth of 2.5 cm, preserving the most superficial layers of the skin and improving its quality.



Enygma

INCORPORATES

· Enygma stimulates the **regeneration and rejuvenation of vulvovaginal tissues**. Due to the **X-Orbital** technology, Enygma can deliver electromagnetic energy at 360° and in a multi-frequency way, providing effective and controlled tissue heating.


· Enygma includes a **system for temperature measurement** with safety-limited parameters, providing a painless therapy without the need of anesthesia.

Indicated for


- Vulvovaginal hyperlaxity.
- Sexual dysfunction.
- Mild prolapse.
- Stress incontinence.
- Reverses symptoms associated with menopause (vaginal dryness).

Benefits


- Non-surgical.
- No anesthesia required.
- Painless.
- Ambulatory.

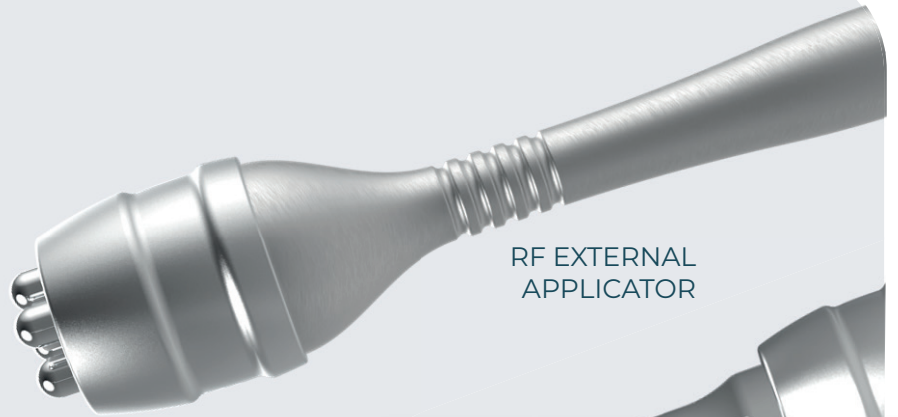

RF EXTERNAL APPLICATOR

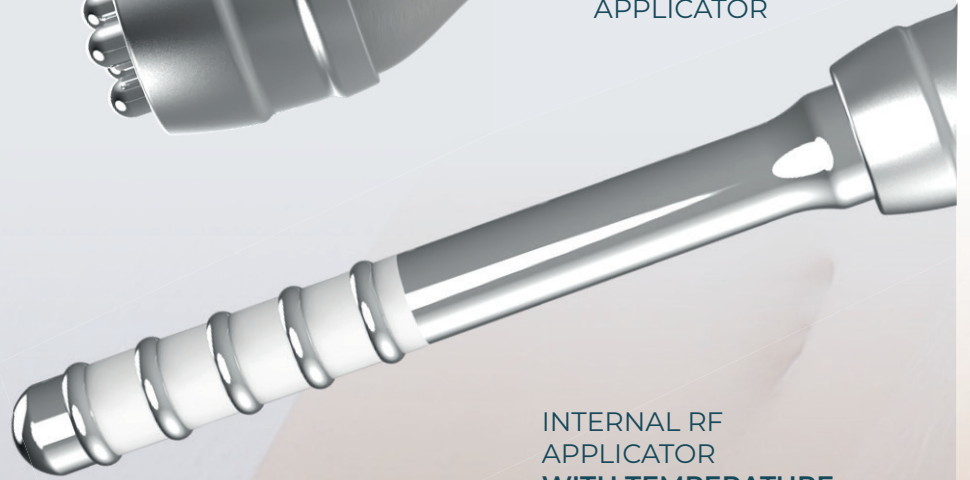

INTERNAL RF APPLICATOR WITH TEMPERATURE SENSOR

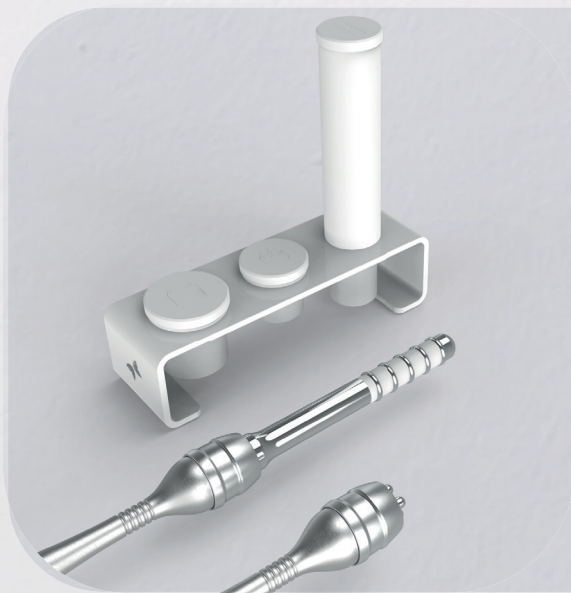

KIT ENYGMA

Add this technology in the following equipment.

- UNYQUE
- ICONYC CRYO + HIMFU



TECHNICAL DATASHEET

COUNTRY OF ORIGIN	ARGENTINA
MANUFACTURER	BODY HEALTH GROUP
PRODUCT NAME	UNYQUE
WIDTH	40 cm
LENGTH	38 cm
HEIGHT	116 cm
WEIGHT	60 kg
CERTIFIED BY	ANMAT
NOMINAL INPUT VOLTAGE	100~240V;50/60Hz
FUSES	T6.3AH, 250V~(5X20mn)
ELECTRONIC SYSTEM	Microprocessor-controlled
ROOM TEMPERATURE	22°C to 28°C microprocessor
HUMIDITY	90% MAX

MAINTENANCE



NO CONSUMABLES



Integral Equipment

Results visible from the first treatment.



Guaranty

1-year full guaranty



User friendly

Easy and efficient PRESET and user interface.



Assembly tutorial



TECHNICAL DATASHEETS

ReFreeze

MAXIMUM VACUUM	200 mBar
FLOW RATE	45 LTS/MIN.
TIME ADJUSTMENT	0 TO 99 MINUTES
ENERGY ADJUSTMENT	1 TO 100%, PITCH 5%
FREQUENCY	500KHz / 800KHz / 1MHz
PART APLIED TO PATIENT	SURGICAL STEEL
TEMPERATURE LIMIT	+5°C / - 5°C / MAX COLD
TYPICAL TEMPERATUR AT MAX	-25°C
TIME ADJUSTMENT	0 TO 99 MINUTES
ENERGY ADJUSTMENT	1 TO 100%, PITCH 5%
CRYO TREATMENT	BODY/FACE
CRYO MODES	BIPOLAR/UNIPOLAR
TREATMENTS W/O CRYO	EYES
MODES W/O CRYO	MULTIPOLAR

MICRODERMOABRASION:

MAXIMUM VACUUM	200 mBar
FLOW RATE	45 LTS/MIN.


CryoFrequency

FREQUENCY	500KHz / 800KHz / 1MHz
PART APLIED TO PATIENT	SURGICAL STEEL
TEMPERATURE LIMIT	+5°C / - 5°C / MAX COLD
TYPICAL TEMPERATURE AT MAX	-25°C
TIME ADJUSTMENT	0 TO 99 MINUTES
ENERGY ADJUSTMENT	1 TO 100%, PITCH 5%
CRYO TREATMENT	BODY/FACE
CRYO MODES	BIPOLAR/UNIPOLAR
TREATMENTS W/O CRYO	EYES
MODES W/O CRYO	MULTIPOLAR

HImFU

OUTPUT	FLOATING ISOLATED BF TYPE (IEC 60601-1)
EMITTER	BARIUM TITANATE
PART APLIED TO PATIENT	SURGICAL STEEL
FREQUENCY OSCILATION	2.9 MHz
OUTPUT ENERGY	45 WATT/CM2
TIME SETTING	0 TO 99 MINUTES
ENERGY ADJUSTMENT	1 TO 100%, PITCH 5%
POWER SETTING SEQUENCE	5 TO 60 SEC ON / 5 SEG OFF
DUTY CYCLE	30/70 40/60 50/50 60/4070/30
FREQUENCY MODULATION	5 / 15 / 30 / 60
COLOR LED	RED
WAVELENGTH LIGHT	624 nm
APERTURE	15°(CADA LED)
LED POWER	320 mW
LIGHT ENERGY	142000 mcd
MAXIMUM POWER	40 mA
LED SEQUENCE	OFF/LOW/HIGH/CONTINUOS




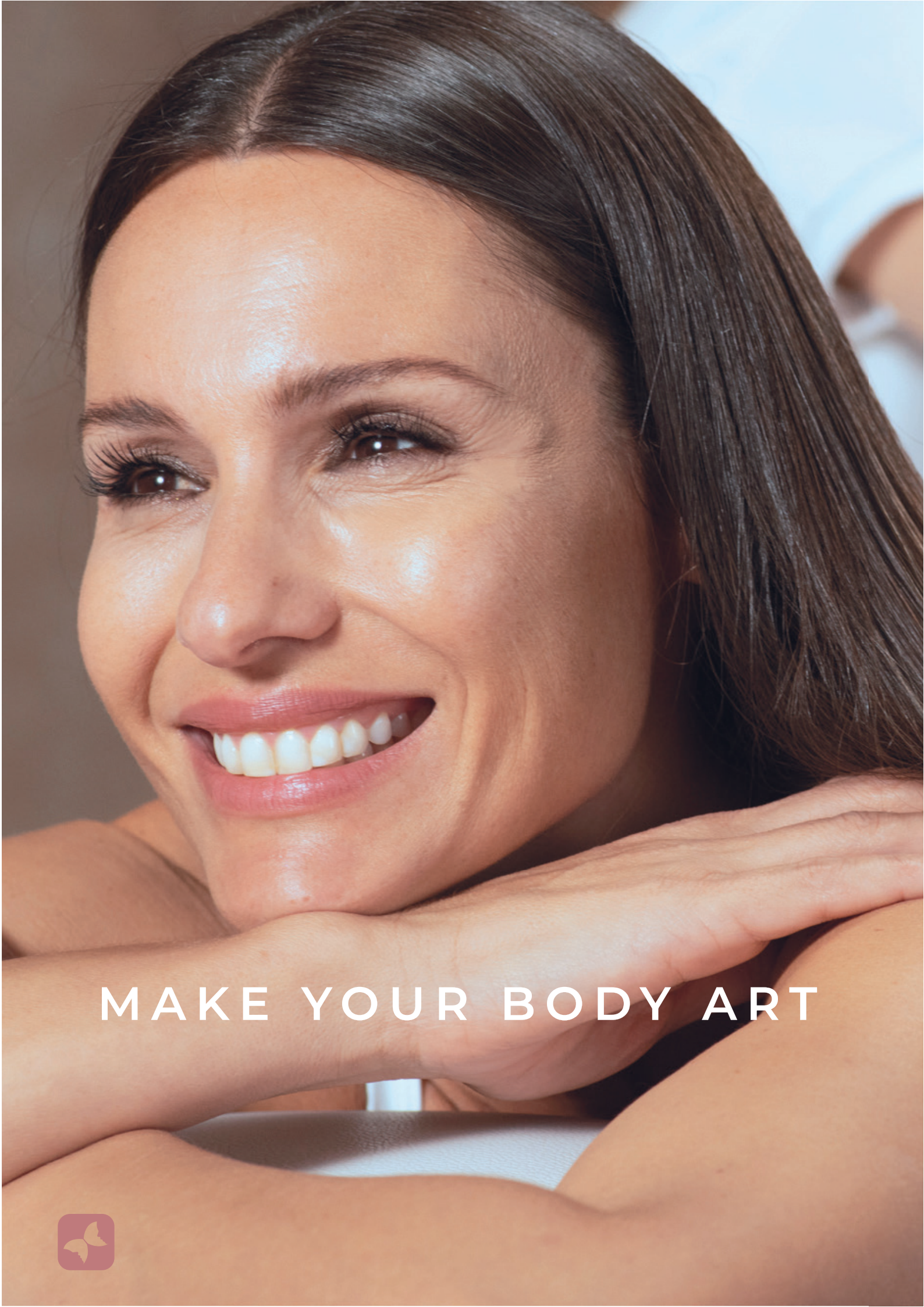

MAKE YOUR BODY ART



AVAILABLE AT 45 COUNTRIES
3 PLANTS IN THE AMERICAN CONTINENT

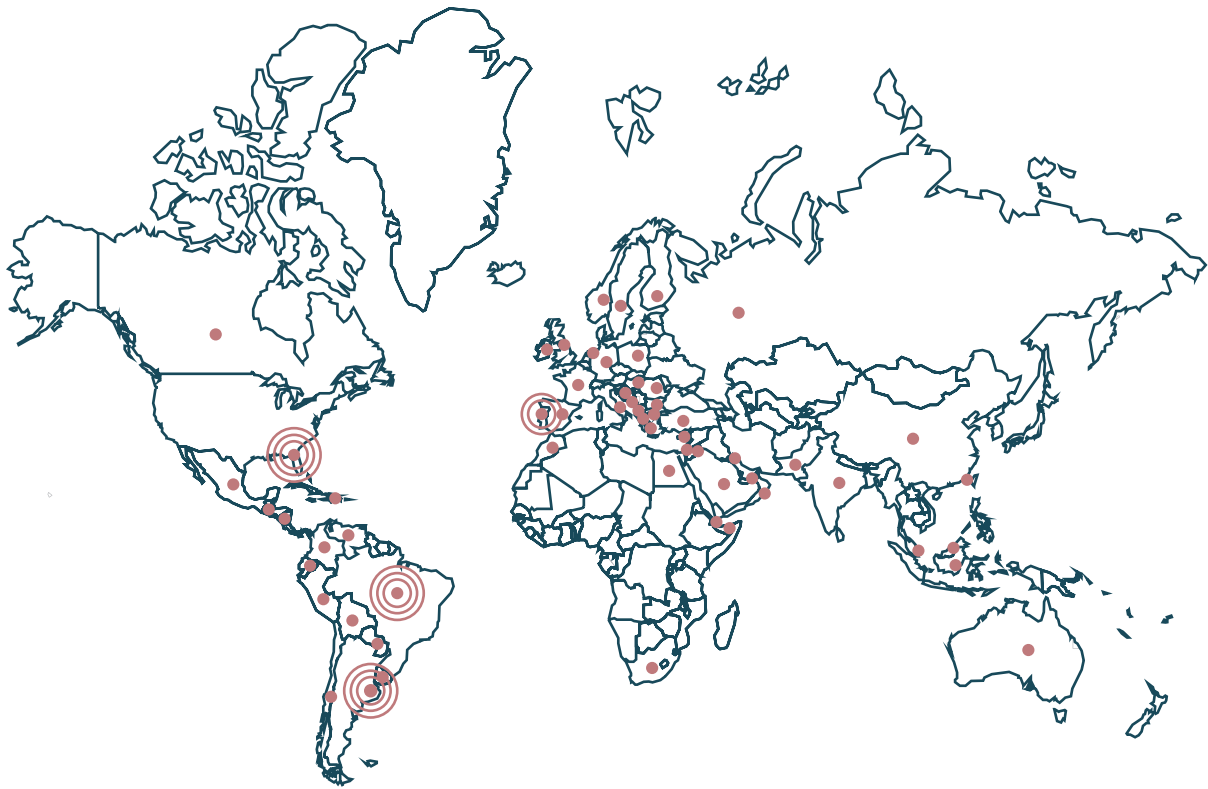

BODY HEALTH ARGENTINA

Belgrano 2638, Torcuato Tigre

Bs. As. - Argentina.

☎ (+54) 11 4727 4747

BODY HEALTH EU

Rua Padre Manuel Valente Pinho Leão 247

4430-684 Oliveira do Douro

Vila Nova de Gaia

Oporto – Portugal

☎ (+35) 19 6145 8471

BODY HEALTH US

10380 W State Rd. 84 Davie Miami

FLORIDA - United State

BODY HEALTH BRASIL

Rua Joaquim Dias Drumond 100 Henrique Nery

Sete Lagoas -MG

☎ (+31) 3775 0709



Equipment with CE marking,
Medical Devices Directive 93/42/EEC.

